

Operating Instructions




Unpack the mat and safely dispose of / recycle the protective plastic bag (**keep away from young children as there is a danger of suffocation!**)



Plug in the mains powered receiver using a suitable mains socket in the room or area of your choice. The receiver has a through adaptor (rated at 13 amps Max) so you won't lose the use of a socket.

If required, set the receiver melody / chime to one of the 36 different types by using the small button on the right-hand side of the receiver above the slide volume switch. If required, adjust the volume of the melody or chime to one of three settings using the volume slide switch on the right-hand side of the receiver.

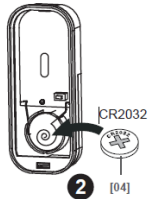

Now remove the battery cover on the mat transmitter unit by inserting a small flat bladed screwdriver into the slot on the end of the transmitter. Insert the coin battery into the mat mounted transmitter unit, positive (+) side facing up, and replace the transmitter battery cover.



With the mat in the same vicinity as the receiver, test the unit by pressing the grey test button on the mat's transmitter unit. A red light will briefly illuminate under the grey test button and the receiver should sound the melody or chime chosen.



Now place the mat in the area to be monitored for foot fall (i.e. besides a bed or under a mat by a door etc). Try to position the mat so that the transmitter unit itself won't be stepped on and subsequently broken. We have provided a non-slip under mat in order to avoid the pressure mat slipping when stepped upon. **Please ensure this is used in order to avoid injury or falls**



Now test the unit again by stepping on the mat. The receiver should again sound the chime or melody confirming that the unit is within range (up to 100m line of site is possible but this will vary according to building layout, construction materials and contents, some experimentation may be required). **The unit is now ready to notify you should someone stand upon the mat.**

Avoid storing the mat with items placed upon it, otherwise the transmitter will be constantly activated and the battery will be depleted when you next come to use it.


Care Instructions




The equipment is for internal use only, keep away from direct heat and sunlight. Use only a damp cloth with a little soap to wipe the pressure mat clean. The mains powered receiver and the battery powered mat transmitter must only be cleaned with a dry cloth and never water or solutions of any kind. The non-slip mat should be cleaned by carefully rinsing the mat under warm running water to remove dust and fibres, this will reactivate the non-slip properties of the mat. To avoid irreparable damage, the pressure mat must not be folded or creased, it must only be used and stored flat. Damage caused to the equipment from improper use / storage are not covered by our warranty. The batteries for the transmitter unit are also excluded from our warranty but we have included a spare transmitter battery for your convenience. Please see the enclosed sheet for details on battery replacement and if necessary, pairing of the transmitter to the receiver unit.

If you should have any questions regards the operation or maintenance of the unit then please email us at info@nursecallequipment.co.uk


Instructions to unpair and then pair the mat transmitter from the mains powered receiver

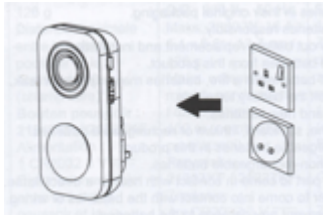

Remove the Mains Powered receiver from the wall socket


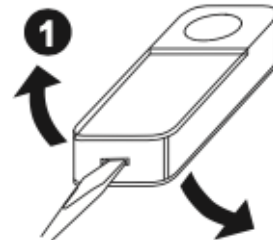


Now remove the battery cover on the mat transmitter by inserting a small flat bladed screwdriver into the slot on the end of the transmitter

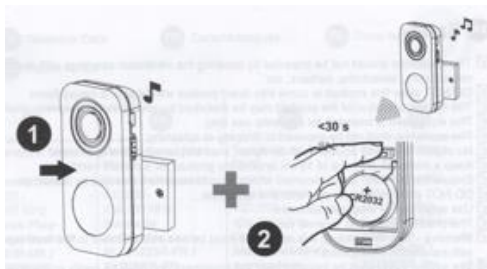

Plug the mains receiver back into the wall socket and within 30 seconds push the small white button mounted internally on the mat transmitter circuit board. The receiver should now chime to indicate it has been unpaired from the mat transmitter. Now replace the battery cover on the mat transmitter. (At this point you should repeat steps 1 and 2 on any other mat and receiver combinations you wish to work separately).

Each mat transmitter and receiver can now be re-paired individually as follows

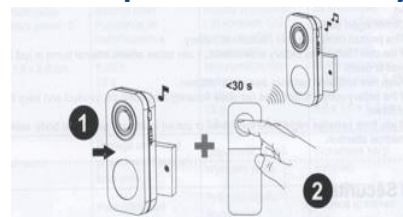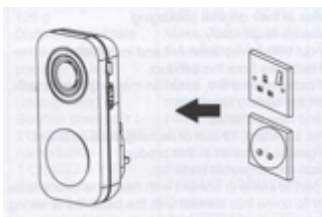

Remove the mains powered receiver from the wall socket and replace again, the receiver will chime, within 30 seconds push the grey test button on the front of the mat transmitter, once again the receiver will chime indicating that this transmitter and receiver have been individually paired. Test the unit by standing on the mat, the receiver should again chime. Please wait for two minutes between pairing any other mat transmitter and receiver combinations, this way all combinations will operate in isolation from each other.