

Operating Instructions




Unpack the mat and safely dispose of / recycle the protective plastic bag
(keep away from young children as there is a danger of suffocation!)



The watch will need to be fully charged by using the USB charger cable included. A red light will illuminate whilst charging. Once fully charged the red charging light will extinguish (2-3 hours) the watch can then be switched on by pressing the power button (bottom button with red circle beside it).



With the mat in the same vicinity as the watch, test the unit by pressing the orange button with the bell symbol on the mat's transmitter unit. A red light will briefly illuminate on the transmitter and the watch should illuminate, sound and vibrate.

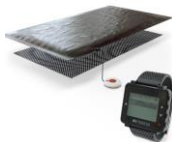

Now place the mat in the area to be monitored for foot fall (i.e. besides a bed or under a mat by a door etc.) Try to position the mat so that the transmitter unit itself won't be stepped on and subsequently broken. We have provided a non-slip under mat to avoid the pressure mat slipping when stepped upon.
Please ensure this is used in order to avoid injury or falls.



Now test the unit again by having someone step on the mat whilst the person with the watch is in the required receiver location. The watch should again illuminate, sound and vibrate confirming that the unit is within range (up to 100m line of site is possible but this will vary according to building layout, construction materials and contents, some experimentation may be required). Briefly pushing the power button twice on the watch (red circle besides it) returns the watch to the normal time display. **The unit is now ready to notify you when someone stands upon the mat.**

Avoid storing the mat with items placed upon it, otherwise the transmitter will be constantly activated and the battery will be depleted when you next come to use it.

Care Instructions



The equipment is for internal use only, keep away from direct heat and sunlight. Use only a damp cloth with a little soap to wipe the pressure mat clean. The mains powered receiver and the battery powered mat transmitter must only be cleaned with a dry cloth and never water or solutions of any kind. The non-slip mat should be cleaned by carefully rinsing the mat under warm running water to remove dust and fibres, this will reactivate the non-slip properties of the mat. To avoid irreparable damage, the pressure mat must not be folded or creased, it must only be used and stored flat. Damage caused to the equipment from improper use / storage are not covered by our warranty. The battery for the mat transmitter unit is also excluded from our warranty, it can be replaced by removing the three small screws from the back of the mat transmitter unit and replace the existing battery with the same type and ensuring that it is fitted the correct way around. Replace the cover and the three screws previously removed and test for correct operation.

If you should have any questions regards the operation or maintenance of the unit then please email us at info@nursecallequipment.co.uk